**LANAFORM**

WELLNESS

# Circulation Mat

ELECTROSTIMULATION  
MAT FOR THE FEET

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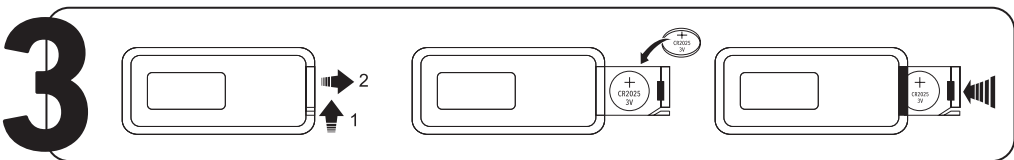
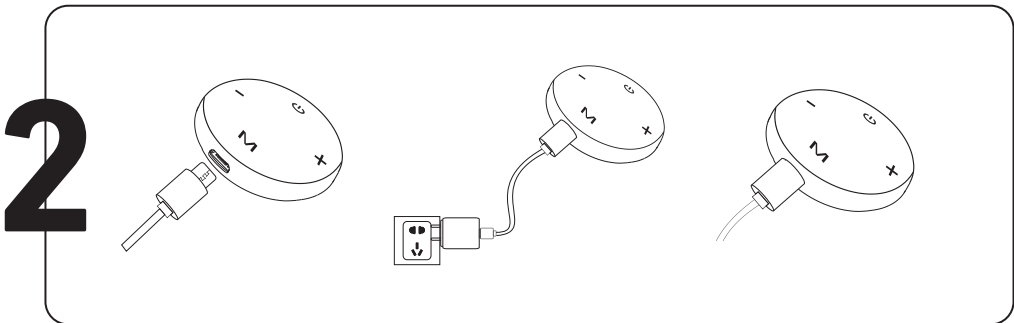
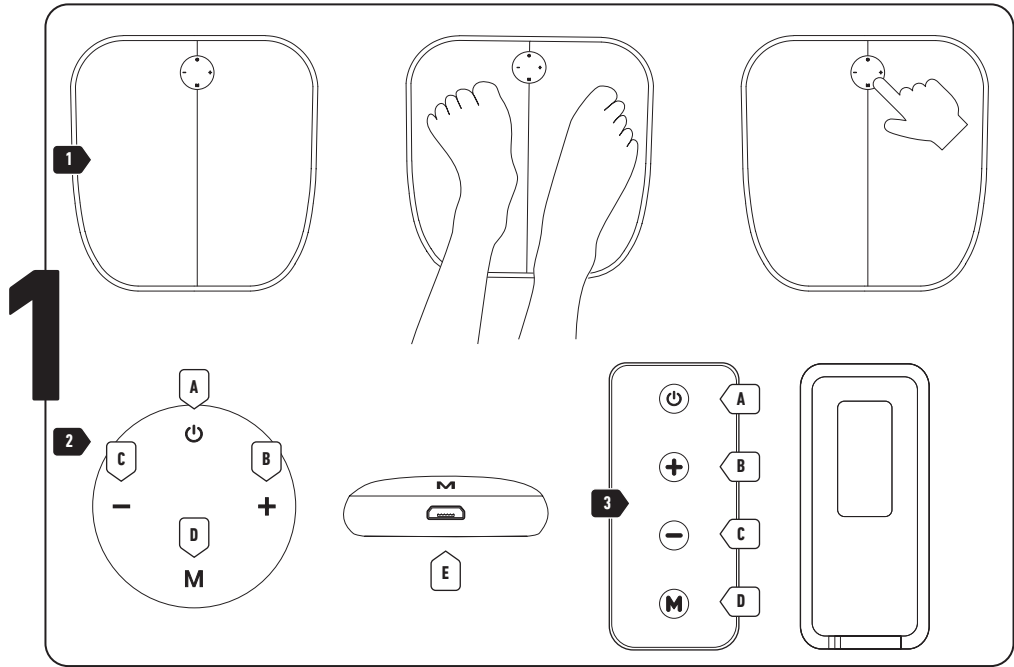
## ENGLISH

### INTRODUCTION

Thank you for choosing the Circulation Mat electrostimulation system by LANAFORM. The Circulation Mat is an EMS electrostimulation mat for the feet for improving blood flow in the lower limbs. Extra fine and light and can be controlled by a remote control, it is very convenient to use anywhere you need it. Rechargeable using a USB cable and foldable, it can be transported very easily in its carrying case.

### Contents of the box

- 1 electro-stimulation device
- 1 main rechargeable unit
- 1 USB power cable
- 1 remote control
- 1 instruction manual

### What is electro-stimulation?

The operating principle of electro-stimulation devices is based on the imitation of the body's pulses transmitted by the skin to the nerve fibres and muscles. The electrical stimuli are safe and painless. The effects of electrostimulation can be quickly felt. We recommend using the Circulation Mat once a day for 25 minutes.

Muscle electrostimulation (EMS) is a recognised technology causing the muscles to contract then relax. Muscle electrostimulation applied to the sole of the feet causes the leg muscles to contract and the veins to compress, which helps send blood in the veins back up to the top of the body, i.e. to the heart. This compression is the natural result of the muscles' action when we walk.

Electrostimulation of the muscles in the lower limbs is used to:

- Tone and relax leg muscles
- Improve local blood flow
- Prevent muscles in the legs weakening and make them mobile in the event of inactivity
- Help soothe symptoms such as pain in the legs and feet, swollen ankles and feet and muscle cramps.

### INSTRUCTIONS FOR USE

**▲ Use of the device is not intended to be a substitute for a medical consultation or treatment. Always consult your doctor therefore beforehand if you are in pain or ill!**

**▲ PLEASE READ ALL THE INSTRUCTIONS BEFORE USING YOUR DEVICE, ESPECIALLY THIS BASIC SAFETY INFORMATION:**

#### **You must consult your doctor before using the device in the following cases:**

- If the symptoms of your pain have not been diagnosed
- Acute illnesses, particularly if high blood pressure, bleeding disorders, tendency towards thromboses and embolisms are suspected or present and in the event of malignant tumours
- All skin diseases
- Chronic, unexplained pain on any part of the body whatsoever
- Diabetes
- Any sensitivity disorders with reduced sensitivity to pain (disorders of the metabolism, for example)
- During other medical treatment
- Pain associated with treatment by stimulation
- Permanent irritation of the skin due to prolonged stimulation using an electrode placed in the same place.

#### **To avoid harmful effects to your health, you must seek medical advice before using the Circulation Mat in the following cases:**

- If you have electrical devices implanted (a pacemaker, for example)
- Metal implants are present
- On people carrying an insulin pump
- If you have a high temperature (> 39°C, for example)
- If you have a known or acute irregular heartbeat and if you have other cardiac excitation and output dysfunctions
- If you have a crisis (epilepsy, for example)
- If you know you are pregnant
- If you have cancer
- After surgical procedures, if stronger muscle contractions might prevent healing

#### **Do not use this equipment in the following cases:**

- On skin with acute or chronic disease (wound or irritation) (in the event of painful and painless irritation, rosacea, itchy skin rashes and hives (allergies, for example), burns, contusions, swellings and open wounds and during the healing process, on scars that are recovering)
- In very humid environments such as the bathroom, in the bath or shower, for example.
- After consuming alcohol
- When a high-frequency surgical device is plugged in at the same time

#### **▲ PLEASE NOTE!**

##### **Use the Circulation Mat only:**

- On men
- For the purpose for which it was designed and in the manner indicated in the instructions for use. It may be dangerous if used incorrectly
- For external use only
- With the accessories originally supplied, which can be ordered optionally (when using other accessories, the warranty will be void)

##### **Precautionary measures:**

- Electro-stimulation is not intended to be used without supervision by people with emotional problems, dementia or a low IQ.
- To avoid injuring the most sensitive skin (extremely rare case), always pull the electrodes gently to remove them.
- Keep the device away from heat sources and do not use it close to (less than 1 metre) shortwave or microwave devices (mobile phones, for example) because this could cause unpleasant spikes in the current.
- Do not expose the device to direct sunlight or high temperatures.
- Protect the device from dust, dirt and humidity. Do not immerse the device in water or other liquids.
- This device is intended for personal use.
- If the device is not working, stop using it immediately in the event of illness or pain.
- To guard against any potential risk, keep out of reach of children.
- Do not use this device at the same time as other devices that send electrical pulses into your body.
- Do not use the device close to highly inflammable materials, gas or explosives.
- Use the device sitting down. If you feel light-headed, stop using the device immediately and raise your legs (for 5 to 10 minutes).
- Make sure that no metal object can come in contact with the mat during simulation. If you wear jewellery or have any piercings (ankle bracelet, for example) in the area of application, remove them before using the device because they might cause localised burns.

##### **Damage:**

- If the device is damaged, do not use it if in doubt and contact your reseller or return it to the address given by the after-sales service.
- Check the device to detect signs of wear and tear or damage. If the device shows signs like this or if it has been used in a non-compliant manner, it must be sent

back to the manufacturer or returned to the reseller before being used again.

- Switch the device off immediately if it is faulty or if it displays malfunctions.
- Do not under any circumstances try to open and/or repair the device yourself! Repairs must be carried out only by the after-sales service or approved resellers. Failure to follow these instructions may void the warranty.
- The manufacturer cannot be held responsible for damage caused by inappropriate or improper use.

#### ▲ Notes on batteries:

- If the liquid of the battery cell comes into contact with the skin or eyes, rinse the area affected with water and consult a doctor.
- Ingestion hazard! Young children might swallow the batteries and suffocate. So keep the batteries out of the reach of children!
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Protect the batteries from excessive heat.
- Explosion hazard! Do not throw batteries onto the fire.
- Do not attempt to re-charge or short-circuit the batteries.
- Remove the batteries from the battery compartment if the device is not to be used for a lengthy period.
- Use only identical or equivalent batteries.
- Always replace all the batteries at the same time.
- Do not use the accumulator!
- Do not dismantle, open or break the batteries.

## DESCRIPTION OF THE DEVICE 1

- 1 Electro-stimulation mat
- 2 Main unit
  - a On/off button
  - b Button to increase the intensity
  - c Button to reduce the intensity
  - d Mode selection button
  - e USB micro port
- 3 Remote control
  - a Off button
  - b Button to increase the intensity
  - c Button to reduce the intensity
  - d Mode selection button

## INSTRUCTIONS FOR USE

#### Ⓞ Notes:

- We recommend using the Circulation Mat once a day for 25 minutes ( automatic timing in each mode).
- The device has 2 modes and 15 levels of intensity.

- We would recommend starting with a weak intensity and increasing the intensity gradually as you adjust it.

## Charging the main battery 2

- 1 Connect the USB cable to the main unit via the USB micro port (Fig. 1, point 3e)
- 2 Connect the USB cable to the mains using an adapter (not included) or to a USB port (for example, on your computer).
- 3 The LED lamp displays a red light when charging. A white light indicates that charging is complete.

## Preparation 1

- Make sure the main unit has been charged (see above).
- Fix the main unit (2) to the mat (1) in the space provided for this purpose.
- Take off your shoes and socks.
- Wash or moisturise your feet. Feet that are dirty or too dry might hinder the feeling of electrical impulses.
- Put each of your bare feet on either side of the mat (1).

## Use 1

- Press the on/off button of the main unit (2a) to switch on the device. (1)

Ⓞ Note: The unit can be switched off by means of the main unit or the remote control, but can only be switched on via the main unit.

- On start-up, the electrostimulation is programmed for mode 1, intensity 1.
- Press the mode selection button, either on the main unit (2d) or on the remote control (3d) to change mode. Le LED light will change colour according to the mode selected: white for mode 1 and green for mode 2.
- Press the button for increasing the intensity either on the main unit (2b) or on the remote control (3b) to increase the intensity. The LED light will flash once.
- Press the button for reducing the intensity either on the main unit (2c) or on the remote control (3c) to reduce the intensity. The LED light will flash once.

Ⓞ Note: if your feet leave the mat, the LED light will flash, the intensity will drop automatically to level 1 and the device will turn off after a few seconds.

## Programmes

Mode 1: adjustment phase lasting twenty seconds or so, then toning programme by intervals.

Mode 2: massage programmes alternating between vibration and kneading.

## Changing the remote control battery 3

- 1 Using your nails, press and pull the battery compartment located on the back of the remote control at the same time, as shown on the diagram.
- 2 Replace the CR2025 battery.
- 3 Push the compartment back inside the remote control.

## MAINTENANCE AND STORAGE

### Cleaning

- Before cleaning, remove the batteries from the device.
- After each use, clean the device with a soft, slightly damp cloth. If it is very dirty, you can use a cloth moistened with a gentle, soapy detergent.
- Make sure that the water does not seep into the device. If this happens, do not use the device again until it is completely dry.
- To clean it, do not use any chemical cleaning agent or abrasive product.

### Storage

- If you do not intend to use the device for a long time, disconnect the main unit from the mat.
- Store the device in a well ventilated area, away from heat and humidity.
- Do not place heavy objects on the device.

## BREAKDOWNS AND SOLUTIONS

### Ⓞ The feeling seems more powerful in one foot than in the other.

Ⓞ This is completely normal. The body reacts differently depending on the areas being treated.

Ⓞ Make sure that your feet are in a good position on the mat.

### Ⓞ When in use, you will experience a burning sensation or stimulation becomes too weak.

Ⓞ Your feet are not sticking properly to the mat.

Ⓞ Apply a few drops of water to the mat or to your feet.

### Ⓞ Stimulation is uncomfortable.

Ⓞ The intensity is too high.

Ⓞ Reduce the intensity.

Ⓞ The mat is damaged.

Ⓞ Contact your supplier's after-sales service.

### Ⓞ The device will no longer turn on.

- Ⓞ The battery is too low.
- Ⓞ Charge the main unit again.
- Ⓞ The remote control battery is flat.
- Ⓞ Use the buttons on the main unit while waiting to replace the battery in the remote control.
- Ⓞ There is a loose connection.
- Ⓞ Check that the main unit is fixed correctly to the mat and that nothing has been damaged.
- Ⓞ **Even after increasing the intensity, stimulation is still too weak or the user cannot feel anything.**

- Ⓞ Is your skin not clean enough?
- Ⓞ Use the device again after cleansing the skin with a damp cloth.
- Ⓞ Are your feet in the correct position?
- Ⓞ Change their position.

If the breakdown persists, do not try to repair or open the device yourself, but return it to the supplier.

## TECHNICAL DATA

- Voltage: battery rechargeable 3.7 V; battery CR2025 3V

- Frequency: 4-35Hz/4-41Hz (with variations)

## ADVICE REGARDING THE DISPOSAL OF WASTE

All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.

When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law.

Prior to disposal, please remove the battery and dispose of it at a collection point so that it can be recycled.

Under no circumstances should spent batteries be disposed of with household waste.

## LIMITED WARRANTY

LANAFORM guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM's control.

LANAFORM may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM's Service Centre. Any attempt to maintain this product by a person other than LANAFORM's Service Centre will render this guarantee void.

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