**LANAFORM**


# Bodyscan Massager


**EN** ENGLISH                      MASSAGING SEAT COVER  
**FR** FRANÇAIS                    COUVRE-SIÈGE MASSANT  
**NL** NEDERLANDS                MASSAGEZETELHOES  
**DE** DEUTSCH                     MASSAGE-SITZAUFLAGE

**PL** POLSKI                            SIEDZISKO DO MASAŻU  
**SL** SLOVENŠČINA                MASAŽNA SEDEŽNA PREVLEKA  
**RO** ROMÂNĂ                        SCAUN PENTRU MASAJ

# 1
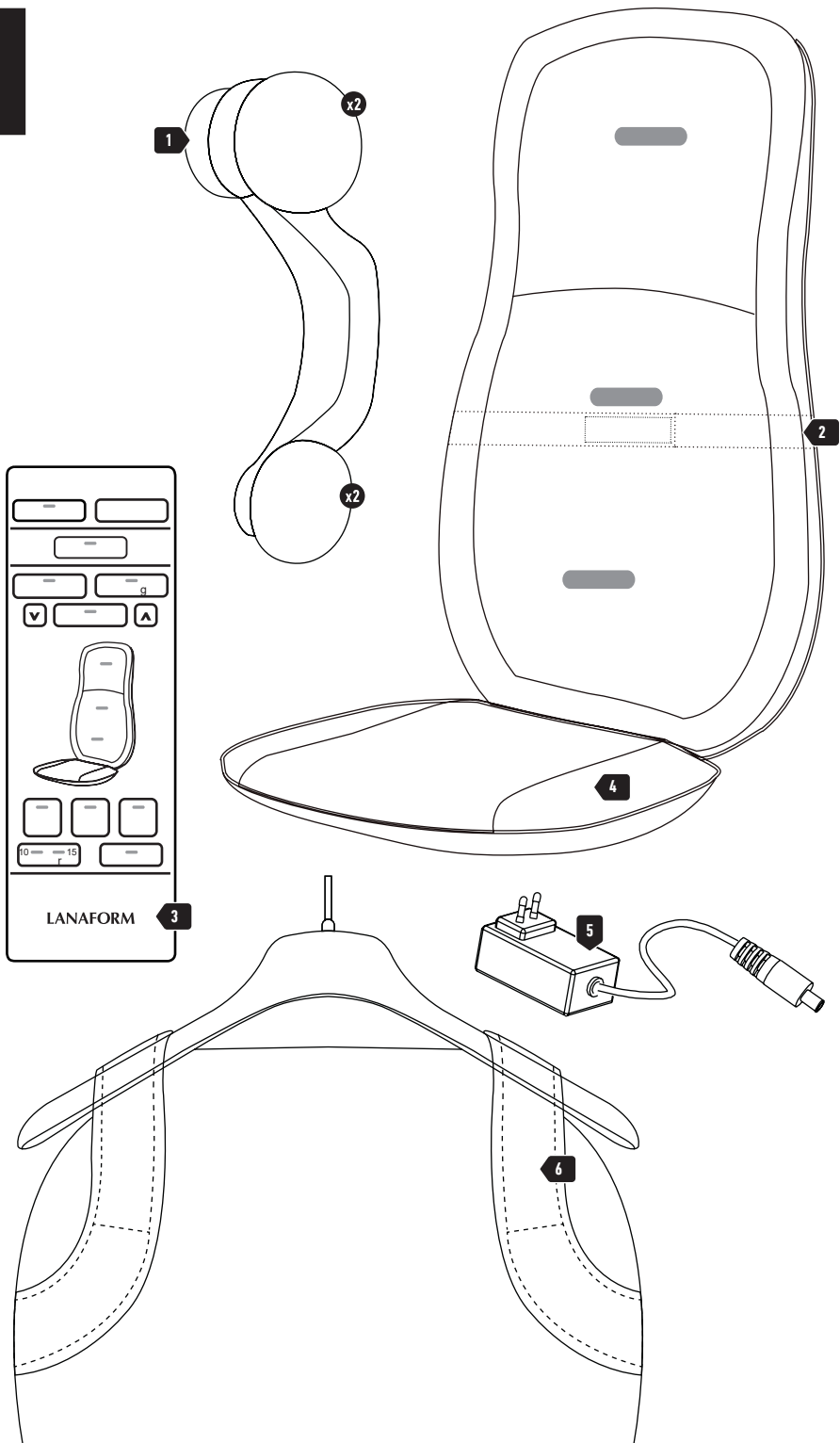
# 2
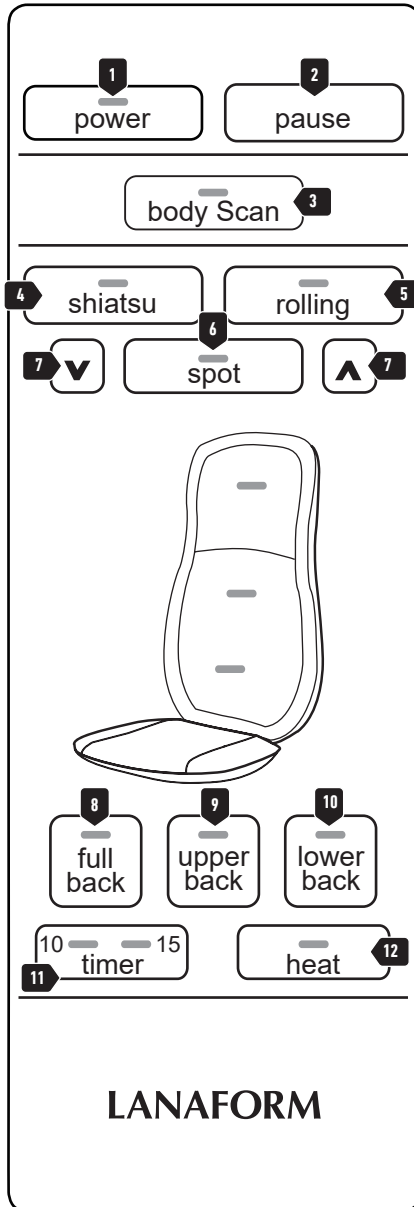
## Contents of the packaging

- 1 Bodyscan Massager seat cover
- 1 remote control
- 1 power cable
- 1 instruction manual

## Precautions for use

**⚠ Please read all the instructions before using your product, in particular these few fundamental safety instructions.**

**📷 The photographs and other product images in these instructions and on the packaging are as accurate as possible but may not provide a perfect representation of the product.**

### **⚠ Warning**

This device can be used by children aged 8 years and over and by people with reduced physical, sensory or mental capacities, or a lack of experience and knowledge, providing they are supervised or have been instructed on the safe use of the device and understand the risks involved. **Children must not be allowed to play with the device. Cleaning and maintenance by the user must not be performed**

by children without supervision. **Children must not play with the device. The surface of the device is hot. People who are sensitive to heat should take care when using this device. Do not use accessories that are not recommended by Lanaform or that are not supplied with this device. If the power cable is damaged, it must be replaced by a similar one available from the supplier or its after-sales department. Do not use this device if the plug is damaged, if the device is not working properly, has been dropped on the floor, if it is damaged or if it has fallen into water. Have the device examined and repaired by the supplier or their after-sales department. Do not carry this device by its power cable and do not use the cable as a handle. Never use drawing pins or other metal fixings with this device. Always unplug the device after use or before cleaning it. If a device falls into water unplug it immediately, before retrieving it. An electrical device must never be**

left plugged in and unattended. **Unplug** it whenever you are not using it. **Keep** the power cable away from hot surfaces. **Do not** use this product before going to sleep. **Massage** acts as a stimulant that may delay the onset of sleep. **Never** allow anything to fall into and never insert any object into any of the openings. **Never** use this device in a room in which aerosols (sprays) are used or in a room in which oxygen is administered. **Do not** use this device under a cover or under a cushion. **Excessive** heat could cause a fire, electrocution or injury to a person. **To** unplug the device, ensure the "POWER" button is switched off before removing the plug from the power socket. **If** you are concerned about your health, consult a doctor before using this device. **Do not** use the device:

- In the event of pathological changes or lesions on the back (e.g. a slipped disc, an open wound);
- While pregnant;
- While sleeping;
- In a vehicle;
- For animals;

- While doing activities during which an unexpected reaction could be dangerous;
- Following the absorption of substances that limit sensory perception (e.g. drugs with an analgesic effect, alcohol).

**If** you experience any pain while using this device, stop using it immediately and consult your doctor. **Do not** use this device for more than 15 minutes at a time (risk of overheating). **Leave** it to cool for at least 15 minutes before re-using. **Never** use this device on any part of the body that is swollen or inflamed or if you have a rash. **Seek** medical advice before using this massage device:

- If you suffer from a serious illness or if you have undergone surgery on your upper body;
- If you have a pacemaker, implant(s) or any other medical device, thrombosis, diabetes;
- If you have any pain of undetermined cause.

EN

5 - ENGLISH

This product is a non-professional massage device designed to relax tired muscles. **Do not use this device as a substitute for medical treatment. It is intended only for home and non-professional use. Never use under a blanket/sheet as this could obstruct the air hole. Excessive heat could occur, causing fire, electric shock, or injury.**

### **⚠ Caution**

Always wear a shirt/clothing during use. **Do not use on bare skin to avoid potential irritation. Do not use on areas of sensitive skin or on areas with poor circulation. Do not expose to extreme temperatures or to humidity. As with any electrical device, this massage device must be used with caution and care to avoid all risk of electrocution. Therefore, the device must:**

- only be used at the mains voltage indicated on the device;
- never be used if the device or its accessories show any signs of damage;
- not be used during a storm.

If there is a fault or malfunction, switch off the device and unplug it immediately. **Do not pull on the power cable or on the device to unplug it from the power supply. Never hold or carry the device by its power cable. Do not squash, fold or twist the cable. Do not push pins or other sharp objects into it. Ensure the massage device, the switch, plug and power cable do not come into contact with water, steam or other liquids. Therefore, the device must only be used:**

- indoors, in rooms without high humidity (e.g. never in the bathroom, in a sauna),
- with dry hands

### **Presentation**

Thank you for choosing Lanaform's Bodyscan Massager. The Bodyscan Massager is a massaging seat cover with 4 massage heads that perfectly imitate shiatsu massage and rolling massage. With its scanning function, you can benefit from a unique massage tailored to your body! Whether it is used for the whole back and neck or just for a specific point, the massaging effect is particularly effective for relieving contracted muscles and toning tired tissues. The device brings you beneficial relaxation that will make you forget your day-to-day stress.

# Description

## The device 1

- 01 Four massage heads
- 02 Fastening straps
- 03 Manual control
- 04 Seat cushion
- 05 Power cable
- 06 Hanger strap

## Remote control 2

- 01 POWER: on/off button and initial reset of the massage mechanism
- 02 PAUSE: stops the device immediately
- 03 BODY SCAN: button to activate the body scan
- 04 SHIATSU: button to select the shiatsu massage mode
- 05 ROLLING: button to select the rolling massage mode
- 06 SPOT: button to select a precise massage area (spot massage)
- 07 UP and DOWN ARROWS: buttons to adjust the height (spot massage)
- 08 FULL BACK: whole back and neck massage
- 09 UPPER BACK: massage for the upper back area
- 10 LOWER BACK: massage for the lower back area
- 11 TIMER: button to select 10 or 15 minutes
- 12 HEAT: button to turn the heating function on/off

## Preparation before use

- 01 Place the cushion upright on a suitable chair with a seat and backrest. Ensure the seat section is placed so that it is in full contact with the seat of the chair and at an adequate height for the back.
- 02 Use the straps to secure the Bodysan Massager.
- 03 Plug the device into the mains.

# Instructions

- 01 Press your back vertically onto the massage device. You will need to sit up straight to ensure accurate palpation of the back. This is made possible by the body scan function, which measures the user's size. At the start, press the massage device carefully against the back. Ensure that you are sitting comfortably on the seat, with your back straight. Ensure you position yourself in the middle so that the massage heads move properly to the left and right of your spine. Check carefully that the massage position feels nice to you and then move your weight gradually towards the massage device. The massage should feel pleasant and relaxing at all times. If the massage is painful or unpleasant, stop it or change the position of the device
- 02 Press the POWER button to turn on the massage device. The LED display on the BODYSCAN button flashes throughout the back palpation. After the palpation, the shiatsu massage programme is activated. The LEDs indicate the massage function currently being used.
- 03 Choose your massage mode:
  - SHIATSU: A combination of shiatsu and rolling massage applied to the back and neck
  - ROLLING: A rolling massage applied to the back and neck.
- 04 Choose the massage area:
  - FULL BACK: the whole back and neck
  - UPPER BACK: the upper back
  - LOWER BACK: the lower back
  - SPOT: using the "up" and "down" arrows, select a specific place on the back for the massage.
- 05 Turn the heat function on or off by pressing HEAT.

- 06** Select the amount of time (10 or 15 minutes). When choosing the massage area, the automatic stop is activated and set to 15 minutes automatically. By pressing the timer button, the automatic stop can be set to 10 or 15 minutes. Warning: once the timer settings have been changed, the countdown starts again from zero. To prevent the device from overheating, select the maximum massage time at the start and do not change this during the massage.
- 07** If you want to stop the massage before the end of the selected time period, press the POWER button. The massage heads will return to their starting position. Do not disconnect the device from the mains during operation. Always wait a few moments after turning off the device to allow it to return to its starting position.

## Cleaning

### ⚠ Maintenance

- 01** Disconnect the cushion from the power source and wait a few minutes before cleaning it.
- 02** Clean it with a soft, slightly damp sponge.
- 03** Never let water or any other liquid come into contact with the cushion.
- 04** Do not immerse it in liquid to clean it.
- 05** Never use abrasive detergents, brushes or thinners to clean it.
- 06** Wait until the device is completely dry before re-using.

## Storage

When not using the device for long periods of time, we recommend that you keep it in the original packaging in a dry place and do not put any items on top of it. You can also hang the device in a wardrobe using a hanger. To do this, the hanger must be inserted into the hanger strap on the side of the device. Make sure you use a stable wooden hanger that can withstand the weight of the device. Avoid contact with sharp edges and other pointed objects that could cut or pierce the fabric surface. Do not hang the cushion up by the manual-control lead or the power cable.

## 🔧 Troubleshooting

🔍 Problem	🕒 Reason	😊 Solution
The massage heads' are moving too slowly.	The massage heads are being subjected to excessive stress.	Reduce the pressure of the massage.
The massage heads are only moving in the upper or lower area.	The massage has only been activated for a certain area.	Press the full, upper back or lower back button to change/extend the area.
The massage heads aren't moving.	The device is not plugged in.	Plug-in the device and turn it on.
	The device is not switched on.	Turn on the massage device using the POWER button. Select the area to be massaged.
	The overheating-protection system has been triggered, the device is damaged.	Contact the retailer or an after-sales department.

## Technical data

- Product reference: LA110318
- Power supply: 12VDC, 2.5A
- Dimensions: 108 x 40 x 10.9cm
- Weight: approx. 9.0 kg

## 🗑️ Waste disposal advice

The packaging is entirely composed of environmentally safe materials that can be dropped off at your local sorting centre for use as secondary materials. Cardboard can be placed in a paper recycling bin. Packaging film should be taken to your local sorting and recycling centre.

Once you are no longer using the device, dispose of it in an environmentally-friendly manner and in accordance with legal guidelines.

## Limited warranty

Lanaform guarantees that this product will be free from any material and workmanship defects for a period of two years from the date it was purchased, other than as specified below.

The Lanaform warranty does not cover damage caused by normal wear and tear of this product. In addition, the warranty for this Lanaform product does not cover any damage caused as a result of misuse or improper use, accidents, attaching any unauthorised accessories, modification of the product or any other condition of any kind that is out of Lanaform's control.

Lanaform shall not be held liable for any type of incidental, consequential or special damage.

All implied warranties of fitness for the product are limited to a period of two years from the original date of purchase, provided that a copy of the proof of purchase can be presented.

Upon receipt, Lanaform will repair or replace your device, as appropriate, and return it to you. A warranty claim can only be made through the Lanaform Service Centre. Any maintenance of this product entrusted to any person other than the Lanaform Service Centre will cancel this warranty.

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