LANAFORM

Cell Roller

EN ENGLISH

SQUEEZE-ROLL MASSAGE

FR FRANÇAIS

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NL NEDERLANDS

PALPATIE-ROLMASSAGE

Presentation

Thank you for choosing the LANAFORM Cell Roller massage device.

The Cell Roller anti-cellulite device makes it possible for you to receive a massage that encourages blood flow, especially in areas affected by cellulite, the orange peel effect and loose connective tissue: the thighs, arms, belly and bottom. Its powerful motor provides a strong squeeze-roll massage.

This complex movement makes it possible to massage the subcutaneous connective tissue, encouraging blood flow and lymphatic circulation. Your skin will become smoother, firmer and cellulite will appear reduced.

Contents

- 01 Unpacking
- 02 Safety precautions
- 03 Product description
- 04 Instructions for use
- 05 Cleaning and maintenance
- 06 Storage
- 07 Advice regarding the disposal of waste
- 08 Limited warranty

Unpacking

Please open the packaging carefully before use. Check whether all the accessories are present and whether any components have been damaged during transport and follow the installation and operating instructions in this manual. In the event of any damage or malfunction, please contact the reseller

or Lanaform directly. In the event of any complaint, you will need the model of the device, serial number, date of purchase, your details and your address.

- 1× Cell Roller device
- 1× Power supply
- 1× this User Manual

Safety precautions

⚠ Please read all the instructions before using your electric hot water bottle, especially this basic safety information.

📷 The photographs and other images of the product in this manual and on the packaging are as accurate as possible, but do not claim to be a perfect reproduction of the product.

The device is intended solely for personal use and not for medical or commercial purposes. It is not a substitute for medical treatment. This appliance can be used by children aged 8 years and over, as well as by those with reduced

EN

4 · ENGLISH

physical, sensory or mental capabilities or lack of experience and knowledge as long as they are supervised or know how to use it completely safely and understand the risks. **C**hildren must not play with the appliance. It must not be cleaned or handled by children without supervision. **I**f the device's electrical power cable is damaged, it must be disposed of. **I**f the cable cannot be detached, the device must be disposed of. **U**nder no circumstances should you open or repair the device; it would no longer be possible to guarantee that the device would work properly. **F**ailure to follow these instructions may invalidate the warranty.

⚠ Warning

The use of a massage device is prohibited

- on animals,

- do not use it if you have any pathological changes or wounds on the part of the body to be massaged (e.g. herniated disc, open wound),
- for massages in the area of the heart if you wear a pacemaker. For other parts of the body, seek medical advice first,
- on skin and swollen or burnt parts of the body presenting with inflammation or injury,
- in the event of open wounds, on capillaries, varicose veins, in cases of acne or rosacea,
- herpes or other skin diseases,
- during pregnancy,
- on the face (eyes), larynx or other particularly sensitive parts of the body,

- do not use the device for a foot massage,
- when sleeping,
- after taking drugs or alcohol (limited sensory perception),
- in a vehicle.

Seek medical advice before using this massage device:

- if you suffer from a serious illness or if you have undergone surgery on your upper body,
- if you have been fitted with a pacemaker, implants or any other ancillary device,
- if you have ever suffered or are currently suffering in particular from blood clots, diabetes or any other pain of unknown origin.

The device is only for use as per this user manual. The manufacturer cannot

EN

7 - ENGLISH

be held liable for damage caused by inappropriate or improper use.

⚠ Electrocutation

As with any electrical device, this massage device must be used cautiously and with care to avoid any risk of electrocution. **Therefore**, only use the device:

- with the power supply provided and at the mains voltage indicated on the power supply,
- under any circumstances if there is visible damage to the device or its accessories,
- during a thunderstorm.

If the equipment is faulty or malfunctions, stop and unplug it immediately.

Do not pull on the device's power cable or to disconnect the power supply.

Do not hold and never carry the device by its power cable. **Keep** cables away

from hot surfaces. **Do** not jam, fold or twist the cable. **Do** not push needles or other sharp objects into it. **Make** sure that the massage device, the switch, the power supply and the cable do not come into contact with water, steam or other liquids. **Therefore**, only use the device only in dry rooms indoors (i.e. never in a bathroom or sauna). **Never** try to pick up a device that has fallen into water. **Unplug** it immediately. **Do** not use the device or its accessories if any damage is visible. **Avoid** shaking the device and do not drop it.

⚠ Repairs

Repairs to electrical appliances should only be carried out by qualified electricians. **Non-compliant** repairs can pose serious risks to the user. **Contact** after-sales service or an authorised dealer for any repairs.

EN

9-ENGLISH

⚠ Fire hazard

Using this appliance incorrectly or failing to comply with these instructions could cause a fire! **T**herefore, never use the massage device:

- under a cover or cushion,
- close to petrol or other highly flammable substances.

⚠ Handling

After each use and before it is cleaned, switch off the device and unplug it. **D**o not insert any item into the openings or into any rotating parts. **M**ake sure that moving parts can always move unimpeded. **K**eep the device away from high temperatures.

EN

Product description 1

- 01 Handle
- 02 Slider switch, 3 positions
 - 0: stop
 - 1: low speed
 - 2: high speed
- 03 Massage rollers
- 04 DC port
- 05 Massage roller connector
- 06 Power supply

Instructions for use

Introduction

Cellulite can cause the formation of fatty deposits on the surface of the skin, located mainly around the thighs, arms and buttocks. When cellulite is widespread, the fat cells increase in volume and thus exert a certain upward pressure, causing this common formation of fatty deposits

on the surface of the skin. To prevent or reduce cellulite, medicine and the world of cosmetics offer many methods of treatment, but yet none of them have yet yielded entirely satisfactory results.

However, massage that encourages blood flow, combined with regular physical exercise and a balanced diet, can lead to an improvement in the appearance of the skin.

The Cell Roller device makes it possible for you to receive a massage that encourages blood flow, especially in areas affected by cellulite orange peel effect and loose connective tissue: the legs, arms and bottom. The powerful motor provides a strong squeeze-roll massage. This complex

EN

11 · ENGLISH

movement makes it possible to massage the subcutaneous connective tissue, encouraging blood flow and lymphatic circulation.

With its squeeze-roll function, the device reproduces the action of a professional massage. The ergonomic design of this device provides for comfortable use at home.

Set-up

- Remove the packaging.
- Check that the appliance, power supply mains plug, and cable are not damaged in any way.
- Insert the connecting plug into the socket on the massage device.
- Position the power cable so that nobody trips over it.

Use

The massage device must not be pressed hard into the area you want to massage. If the massage is painful or unpleasant, stop or change the position of the device.


Do not insert or force any part of your body, especially not your fingers, between the rotating massage rollers.

Turn on the massage device using the slider switch. The massage rollers will start turning. You can choose between low- or high-speed massage. We recommend starting with low-speed. To turn the device off, set the slider switch to the 0 position.

Use on the belly and waist

Sit or stand comfortably with your back straight. Put the Cell Roller over the abdomen and gently tense your belly. Perform slow movements over your entire belly, letting the device work on each area. Starting from the middle of your belly, you can also work in a straight line towards either side. To work the tops of the hips, make forward-and-back movements on each side as shown in the diagram.

Repeat one of these movements several times. After several uses, you can also choose the higher intensity for improved results in less time, depending on your skin and your feeling when using the Cell Roller.

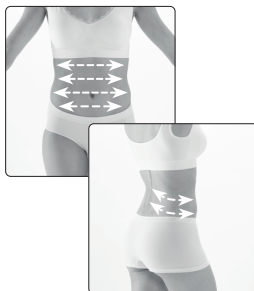⚠ Warning


If you experience any discomfort, stop using the Cell Roller immediately and consult your doctor.

Use on the legs

For an effective massage of the legs, bend over slightly and place the Cell Roller just above the knees. Perform slow movements on the front and back of your thighs until you reach the groin and hips. Alternate between one leg, then the other. The massage should not last longer than 5 minutes per leg.


EN

13 · ENGLISH

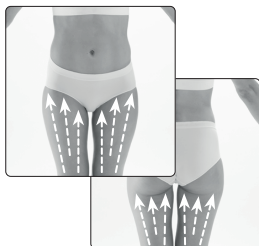Use on 'saddlebags'

More specifically, to work on the sides of your thighs, perform small horizontal and vertical movements. Do not go over 5 minutes per area.

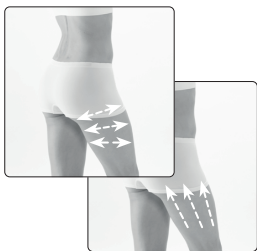Use on the buttocks

To tone your buttocks, start using the Cell Roller at the base of the buttocks and work up towards the waist. Let the device pass over the entire area you want to work on. You can also make horizontal outward movements from the middle. The recommended duration of use is 5 minutes on each buttock.

Use on the arms

To tone your arms, make up-down movements on

the front and back of your arms.

Cleaning and maintenance

⚠ Warning

- Before cleaning, disconnect the cable from the power supply and unplug the device.
- Clean the device using only the method shown. Under no circumstances should liquid

enter the device or its accessories.

- You can remove the massage rollers by pulling them sharply. The rollers can be cleaned with a soft brush, running water and a little detergent.
- The device can be cleaned with a soft damp cloth.
- Wait until the device is completely dry before re-using it.

Storage

- If you are not going to use the device for a long time, we recommend storing it in its original packaging in a dry place. Do not put anything on top of it.

EN

15 - ENGLISH

Advice regarding the disposal of waste

All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.

When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law.

Limited warranty

LANAFORM guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM's control.

LANAFORM may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM's Service Centre. Any attempt to maintain this product by a person other than LANAFORM's Service Centre will render this guarantee void.

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