LANAFORM®

Visage+

Anti-wrinkle device

Appareil anti-rides

Anti-rimpelapparaat

USER'S GUIDE - Visage+


Thank you for choosing the Lanaform «Visage+», a new-generation facial care product for home use. The Lanaform«Visage+» offers a unique and powerful regenerating treatment that makes your skin look younger.

In addition, the «Visage+» helps strengthen tissue fibres thanks to its gentle skin application. It also stimulates blood circulation, helps feed and oxygenate skin tissues, and hydrates the skin, giving a new glow to your face.

The appliance is ergonomic and operates efficiently with 2 AA batteries (not supplied).

APPLIANCE DESCRIPTION:

1. Filter
2. Battery cover
3. Treatment head
4. On/off button

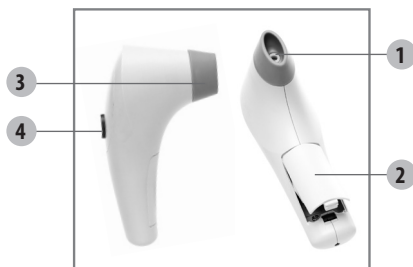1. SAFETY INSTRUCTIONS

Please read these instructions before using the appliance. Please also follow the instructions carefully when using the appliance.

To avoid the risk of burns, electrocution, fire or injury:

1. Only use the appliance in accordance with the instructions for use contained in this manual.
2. The appliance should not be used by persons (including children) whose physical, sensorial or mental capacities are reduced, or who lack experience or knowledge, unless a person responsible for their safety can supervise them or give them preliminary instructions for product use. Children should be supervised to ensure that they do not play with the appliance.
3. Do not use accessories that are not recommended by LANAFORM or that are not supplied with the appliance.
4. If it does not function correctly, falls on the floor, becomes damaged or falls in water, return

- the appliance to the manufacturer or its after-sales service for inspection and repair.
5. Do not use the product before going to sleep. The massage has a stimulating effect that can keep you awake for a period of time.
 6. Never use the appliance in a room in which aerosol products (sprays) are used, or oxygen is administered.
 7. Ensure that the product is not damaged before using it.
 8. The «Visage+» is a non-professional massage appliance. Do not use the appliance in place of medical treatment.
 9. Keep the appliance moving when using it.
 10. Regular cleaning of the filter is necessary to maintain suction efficiency.
 11. The appliance is not waterproof and should not be used in the shower.
 12. Use the «Visage+» on clean and dry skin.
 13. Remove the batteries from the device if you are not using it for a long period.
 14. Never operate the appliance if it has fallen into water.
 15. Do not operate the appliance under a cover or cushion. The appliance might overheat and cause a fire, electrocution or injury.
 16. Do not insert any objects into the appliance's openings.
 17. If you are concerned about your health, consult your doctor before using the appliance.
 18. If you feel any pain when using the appliance, stop using it immediately and consult your doctor.
 19. Do not use the product directly on an inflammation, swellings or skin rashes.
 20. Long periods of use may cause overheating and shorten the appliance's product lifetime. We therefore recommend using the appliance for 10-minute periods, and that you allow it to cool between sessions.

2. PREPARING THE APPLIANCE FOR USE:

Inserting the batteries:

- Make sure that your hands are dry when handling the batteries.
- Place 2 AA batteries (not supplied) in the appliance according to the indicated pole positions.

3. USING THE VISAGE+:


To obtain improved results, remove all make-up, and clean and dry your face.

Use the «Visage+» daily (for approx. 10 minutes) in the morning or evening.


For best results, we recommend you apply an anti-ageing cream or lotion on the face and neck before treatment.

Press the ON / OFF button to activate the appliance.

Move the appliance gently from the centre to the outer edges of the face according to the following instructions:

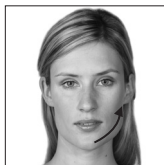
A. Application from the chin to the lower ear area.

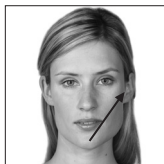
B. Application from the chin to the upper ear area.


C. Application from the corner of the mouth to the sideburn area.

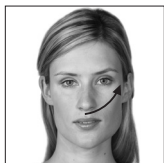
D. Application from under the nose to the sideburn area.

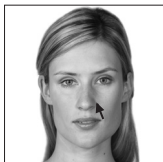
E. Application from the corner of the mouth to the side of the nose.

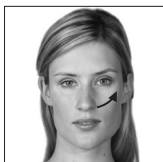F. Application from the side of the nose to the cheekbone.

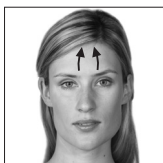G. Application from the bottom to the top of the forehead.

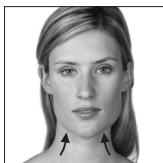H. Application from the bottom to the top of the neck.

I. Application from the bottom to the top of the eye.

J. Application from the middle of the nose to the edge of the eye.

4. CLEANING AND STORAGE:

- Clean the rubber tip after each usage. Use a dry cloth to wipe off any residues of cream.
- Always use a damp cloth to clean the appliance; do not use alcohol, acetone, abrasive products etc.

Cleaning the filter:

Step 1: Remove the rubber tip.

Step 2: Remove the filter protection (small white cover) with tweezers, and turn the appliance over to obtain the small metallic filter in your hand. Rinse and clean both.

Step 3: Dry with a towel and place them back in the appliance.

1 Filter protection; 2 Filter

Note:

Be careful not to mislay the filter and its cover while cleaning them because they are small parts.

Retain the product storage box for future usage.

5. MAINTENANCE AND REPLACEMENT:

If the suction force decreases and is not adequate for use, replace the batteries. In addition, weak suction force can be caused by a clogged filter. Clean the filter before using the appliance again.

RECOMMENDATIONS CONCERNING ELIMINATION OF WASTE

The packaging consists exclusively of materials that are harmless for the environment and can be left with your local waste sorting centre to be used as recycled materials. The cardboard can be put in the paper collection container. The packaging films should be given to local waste sorting centre.

When you no longer intend to use the appliance, discard it in an environment-friendly way, in keeping with local regulations. Please remove the battery first, and turn it into a battery recycling station. Warning! Used batteries should never be discarded with household waste

LIMITED WARRANTY

LANAFORM guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM's control.

LANAFORM may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM's Service Centre®. Any attempt to maintain this product by a person other than LANAFORM's Service Centre will render this guarantee void.

The photographs and other images of the product in this manual and on the packaging are as accurate as possible, but do not claim to be a perfect reproduction of the product.

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