
**LANAFORM**

BEAUTY

Limpia
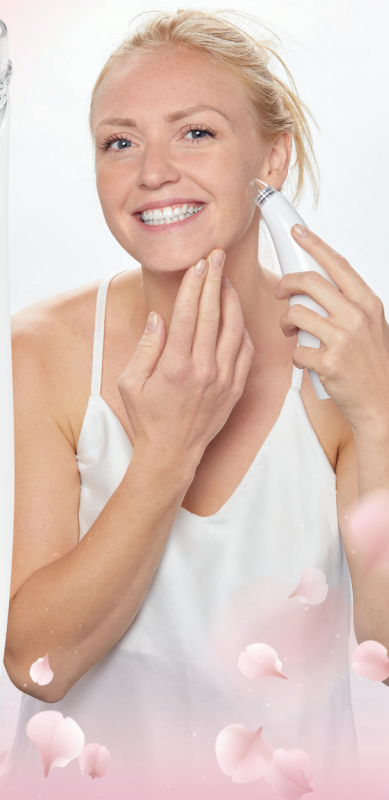
# 1

2

3

## ENGLISH

### INTRODUCTION

Effectively remove pimples, blackheads and impurities with Limpia. This blackhead suction device deeply cleanses the skin, promoting healthier-looking, rejuvenated facial skin with the suction mechanism. It provides individual care and cleansing suited to your needs.

### PRECAUTIONARY STATEMENT

**▲ Please read all the instructions before using your electric hot water bottle, especially this basic safety information.**

 The photographs and other images of the product in this manual and on the packaging are as accurate as possible, but do not claim to be a perfect reproduction of the product.

### ▲ Danger

This device is not intended for use by people, including children, with reduced physical, sensory or mental capabilities, or people with a lack of experience or knowledge, unless someone responsible for their safety can supervise or instruct them beforehand on how to use the device. Children should be supervised to make sure they do not play with the device.

Check that your mains voltage matches that of your device.

Do not connect or disconnect the plug or electrical transformer of the device with wet hands.

If the power cable is damaged, it must be replaced by a similar cable available from the supplier or the supplier's after-sales department.

Electrical appliances should never be left unattended when connected. Unplug it when you are not using it.

Never use this device in a room in which aerosols (sprays) are used or in a room where oxygen is administered.

### ▲ Warning

Do not immerse the device in water or other liquids, and consult the cleaning instructions in this manual.

Keep the device away from heat sources such as stoves, radiators, etc.

Do not leave the device in direct sunlight.

Do not charge the device while in use or if it is wet.

Unplug the device before cleaning it.

Keep the power cable away from hot surfaces.

Never use this device in a room in which the temperature is higher than 40°C.

The lithium-ion rechargeable battery must not be disman-

tled, heated, burnt or exposed to water.

Do not charge the battery prior to use!

When using the device, water or, for example, cream residues from the skin must not penetrate the device because this might cause it to malfunction.

## **▲ Precautions**

Use this device only according to the instructions given in this manual.

This device is exclusively intended to treat the adult human face (except for the mobile eye and lip area).

Do not use the device on animals!

All jewellery (necklaces, earrings, etc.) and contact lenses must be removed prior to use.

Do not use this device to support or replace medical applications..

Never use this device on a body part that is swollen, chapped or inflamed or if you have a skin rash.

The device must not be used by patients with allergies, skin diseases (e.g. herpes, vascular disorders or skin neoplasms), dermatitis, swellings, eczema and other serious (skin) injuries.

Do not use the device on areas of the skin that have been modified by cosmetic surgery.

Do not use the device immediately after exposure to the sun or on sunburn.

This device must not be used on sensitive areas of the skin (e.g. eyes, mouth, ears, neck, etc.) - there is a risk you may be injured!

If you feel any pain at all when using this device, stop using it immediately and consult your doctor.

If you are concerned about your health or are having medical treatment, consult a doctor before using this device.

Do not use this device without medical advice in the following cases: if you are pregnant or breastfeeding, if you have diabetes, if your skin has any tattoos or beauty spots, if you have skin cancer, if you have any unusual skin colouring, if you are taking any medicine.

It may be dangerous if used inappropriately!

This device is designed for domestic use only.

## **COMPONENTS OF THE DEVICE 1**

- 1 Position for interchangeable heads
- 2 Stop/Start button
- 3 Suction power control button
- 4 Light intensity level indicators (3 levels) and battery charge level
- 5 USB charging connector
- 6 USB charger cable
- 7 Normal head

- 8 Large head
- 9 Small head
- 10 Oval head
- 11 Microdermabrasion head
- 12 Replacement filters
- 13 Replacement sealing rings

## INSTRUCTIONS FOR USE

### Skin test

Carry out a skin test 24 hours before using the device. Test the device at the lowest intensity level over small areas of the neck. If the areas tested are reddish or sore after 24 hours, do not use the device.

### Use 2 & 3

 Before starting to use the device, make sure that your facial skin is clean, dry and free of any residue such as powder, cream of makeup.

Your Limpia is ready for use. We recommend using it once a week for combination and greasy skin and once a fortnight for normal and sensitive skin.

- 1 Having cleansed your skin, open your pores with a hot compress or facial sauna.
- 2 Select the head you want and place it on the suction head.
- 3 Press the STOP/START button to start the Limpia and test the suction intensity on the back of your hand. If necessary, adjust the intensity by pressing the intensity control button before using it on your face.

 The device always starts up on the least powerful suction setting.

- 4 Follow the instructions below depending on the area to be treated ( 3).

 The application pressure must be pleasant and not become painful. A slight redness may be visible after treatment. .

**Forehead:** Move the device slowly and uniformly from the centre of the forehead towards the outside. Repeat this movement until this part of your forehead has been treated in its entirety.

**Nose:** Start from the top of your nose and move downwards. Repeat this movement on the right and left sides of your nose.

**Cheeks:** Approximately 2cm away from your eyes, move the device towards the outside of your face, working gently down your cheeks.

**Chin:** Put the device just below your lower lip and work down to the tip of your chin.

### Caution

- Do not use the device in the same place for more than 2 seconds. Carry on moving it over your face in linear movements.
- Make sure you do not move the head more than once over the same area of your face.
- Limit the length of time you use it to 5 minutes.

- You face may have slightly red patches just after treatment. This is normal and the redness will disappear by itself after 30 minutes at the most.
  - Do not use the suction on sensitive skin. Do not use on the fragile skin underneath your eyes.
  - Avoid sweaty activities such as the sauna, solarium or exposure to the sun for one day before and after treatment.
- 5 When you have finished, carefully cleanse your face.
  - 6 Apply a moisturiser or a mask to calm your skin.

## The heads 1

- Normal head (7): This head provides the ideal suction surface for large areas such as the cheeks or forehead. Particularly suitable for normal to combination skin.
- Large head (8): This head is intended for greasier skin, and its even wider suction area makes it possible to remove excess sebum. It can be used over the entire face except on the skin around the eyes.
- Small head (9): This head is suitable for treating areas that are more difficult to access (nose, chin). It is also suitable for sensitive skin because its suction is less powerful than the normal head (7).
- Oval head (10): This accessory is particularly suitable for smoothing wrinkles, promoting elasticity of the skin and blood circulation in the face.

- Microdermabrasion head (11): This head is suitable for peeling facial skin. It exfoliates the skin and removes dead skin. Should not be used too often.

## CHARGING

You can charge your device using the USB cable provided. Connect the USB charger cable to a suitable mains adapter and the device. Then plug the mains adapter into a suitable power socket. You can also charge the device using the USB charger cable through the USB socket on your computer.

The red light will flash while the device is charging. The light will go green once the device is charged.

## MAINTENANCE AND STORAGE

### Maintenance

- Replace the filter each time it is used.
- Wipe the unit with a soft, dry cloth only.
- Disinfect the heads with a conventional disinfectant after each use.
- Never use solvents to clean the device.

### Storage

Having cleaned and dried the device, store it in a dry, cool place out of the reach of children.

## TECHNICAL SPECIFICATIONS

- Reference: LA131311
- 5 interchangeable heads
- 3 levels of intensity
- 10 replacement filters
- 5 replacement sealing rings
- Rechargeable battery: 450 mAh
- Battery life: 65 mins
- Charging time: 60 mins
- Power: 3W
- Voltage: 5V
- Dimensions: 19 x 4.2 cm
- Weight: 290g

## ADVICE REGARDING THE DISPOSAL OF WASTE

All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.

When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law.

Prior to disposal, please remove the battery and dispose of it at a collection point so that it can be recycled.

Under no circumstances should spent batteries be disposed of with household waste.

## LIMITED WARRANTY

LANAFORM guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM's control.

LANAFORM may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM's Service Centre. Any attempt to maintain this product by a person other than LANAFORM's Service Centre will render this guarantee void.

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