LANAFORM

HEALTH

Lumino Plus

LIGHT-THERAPY DEVICE

SAME BENEFITS
AS NATURAL LIGHT!

14 000 LUX
INTENSITY

100 %
UV-FREE

3 × 36 W
TUBES






LAMPE DE
LUMINOTHÉRAPIE

LICHTTHERAPIE-
LAMP

LICHTTHERAPIE-
LAMPE

LÁMPARA DE
LUMINOTERAPIA

LAMPADA DA
LUMINOTERAPIA

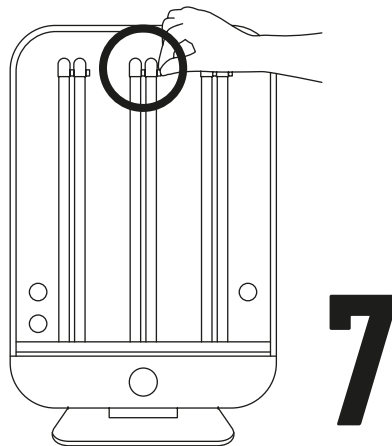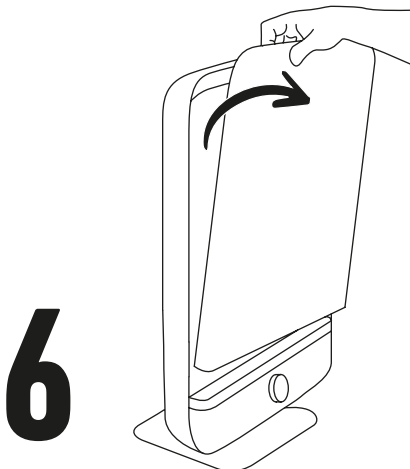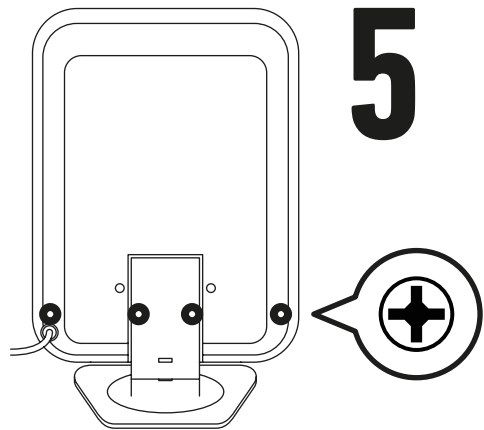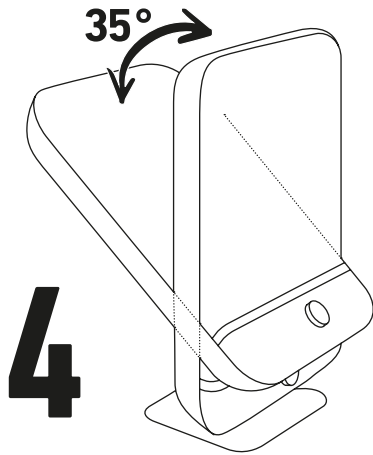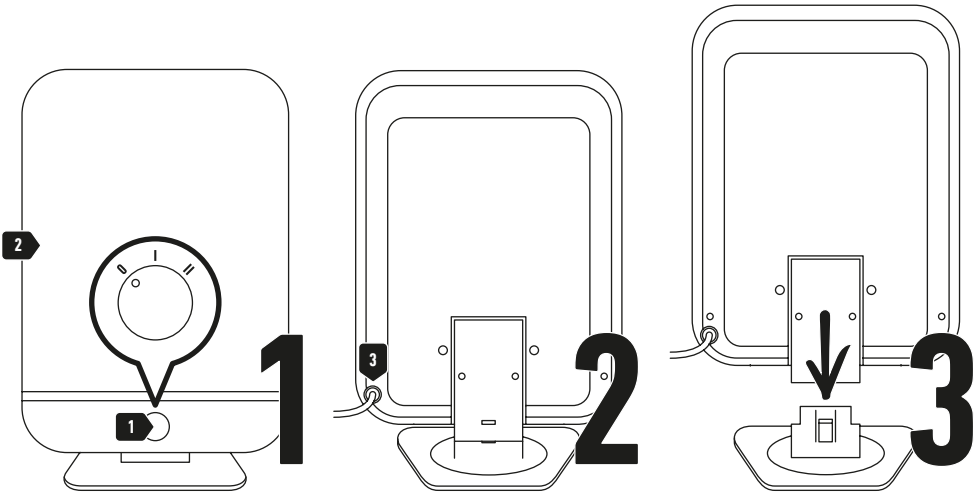ENGLISH

Thank you for purchasing a LANAFORM Lumino Plus. This appliance will help you combat fatigue, apathy and the lack of energy that we often have in autumn and in winter when it is grey outside and the days are getting shorter. These symptoms express what is termed the winter blues.

Lumino Plus is based on the use of a very bright light which simulates the natural rays of the sun. The light produced by Lumino Plus naturally re-establishes the light signals that are essential for feeling good and full of energy.

Thanks to its technology and a higher diffusion of light (about 14000 LUX at 16 cm)*, Lumino Plus will quickly give you your energy and vitality back.

You can use Lumino Plus while doing other activities, such as reading, writing, using the telephone, and even professional activities.

* More than 14000 LUX at a distance of 16 cm when the 3 tubes are lit. In this respect, please refer to the instructions for use.

ABOUT EXPOSURE TO LIGHT

Light and its intensity have a major influence on the body's vital processes. During the autumn and winter seasons, there is less sunshine. The conventional artificial lighting in buildings is not able to offset this deficit in natural light. The consequences of this decrease in luminosity can result in a lack of drive and energy, insomnia or melancholy, problems in concentrating and difficulties getting up in the morning. We refer to all of these symptoms as a whole using the expression "seasonal affective disorder or winter blues".

Many health professionals agree that exposure to luminosity is beneficial in treating the winter blues. By luminosity, they mean daylight with a certain intensity. In winter, the lighting intensity may not exceed 2500 LUX*, while in summer it can reach 100000 LUX on a nice summer day.

LANAFORM has therefore developed this special lighting device that allows you to benefit from the amount of light you need to chase the winter blues away.

Exposure to light consists in using light with a high lux* coefficient in order to simulate natural light from the sun. With regular exposure to light with a high lux coefficient over a given period of time, the symptoms of the winter blues can be greatly reduced, and even disappear.

* The lux is the unit of measurement for the luminance of light.

IMPORTANT SAFETY MEASURES

▲ Please read all the instructions before using your product, particularly these basic safety instructions.

📷 The photographs and other images of the product in this manual and on the packaging are as accurate as possible, but do not claim to be a perfect reproduction of the product.

Only use this appliance for the purposes described in this manual.

This appliance is not intended for use by persons (including children) whose physical, sensorial or mental capabilities are impaired, or by persons lacking experience or knowledge, unless they are supervised or instructed in the use of the appliance first by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.

Before connecting the appliance, check that the voltage indicated on the appliance corresponds to the mains supply voltage.

Please do not use accessories which are not recommended by LANAFORM or which are not supplied with this appliance.

It is prohibited to modify this appliance in any way.

If the power supply lead is damaged, it must be replaced by a similar one available from the supplier or its after-sales service.

Do not use this appliance if the plug is damaged, if it is not working correctly, if it has been dropped or if it is damaged or if it has fallen into water. If this is the case, have the appliance inspected and repaired by the supplier or its after-sales department.

Do not transport this appliance by carrying it by its power lead and do not use this lead as a handle.

Always disconnect the appliance from the mains after use or before cleaning it.

An electrical appliance should never be left plugged in and unsupervised. Disconnect it from the mains if you are not using it.

Keep the power lead away from hot surfaces.

Do not use this product before going to bed. Exposure to light has a stimulating effect which may inhibit sleep.

Never let anything fall into and never insert anything into one of the openings.

Do not use this appliance in a humid environment (in a bathroom, near a shower, etc.).

Never use this appliance in a room where aerosol products (sprays) are used or in a room where oxygen is administered.

Do not use this appliance under a cover or under a cushion. Excessive heat may cause a fire, electrocution of the individual or other injuries.

This appliance is intended exclusively for family use. The use of Lumino Plus is therefore not equivalent to a medical treatment.

Do not use the appliance in a room which is not lit by another source of light, as this can lead to eyestrain. Use this appliance as a supplement for normal lighting (daylight or electric lighting).

If you feel any pain while using this appliance, stop using it immediately and consult your doctor.

We recommend that you do not look at the lamp during your session. Pain in the eyes may occur during the first sessions; this will disappear afterwards.

Consult your doctor before using Lumino Plus if:

- you are suffering or have suffered from serious depression;
- you suffer from hypersensitivity to light;
- you are afflicted with an eye disorder;
- you are taking certain medications (e.g. certain antidepressants, psychotropic drugs, tablets for malaria) which could have an effect contradicting the use of Lumino Plus;
- you have high blood pressure.

Exposure time depends in particular on the distance between you and the Lumino Plus. Refer to the chapter on "Instructions for use" for the recommended distances and exposure times.

COMPOSITION OF THE PRODUCT [1 & 2]

- 1 **ON/OFF** button
 - I Stop
 - II 2 fluorescent economic 36W PLL tubes operating.
 - III 3 fluorescent economic 36W PLL tubes operating.
- 2 Reflector cover
- 3 Electric cable

INSTRUCTIONS FOR USE

It is recommended to use Lumino Plus preferably between 6 am and 8 pm.

Installation

- 1 Insert the base of the Lumino Plus into the main part as shown in the drawing [3]. Make sure that the lamp

- is properly connected to its base before placing it on a solid surface.
- Place your Lumino Plus on a solid surface, such as a table, and adjust the orientation of the Lumino Plus. It can be tilted to 35°. [■ 4]
 - Connect the power lead to the mains.

Use

You can choose from among two type of light intensity:

- Position ①: 2×36W PLL tubes for an intensity of about 7200 LUX at a distance of 20 cm.
- Position ②: 3×36W PLL tubes for an intensity of about 12000 LUX at a distance of 20 cm.

The recommended exposure time to Lumino Plus depends on the distance in relation to the lamp as well as the light intensity selected. However, it is generally recommended to use the appliance for 2 hours at a light intensity of 2500 LUX. According to the table below, you can easily calculate the recommended exposure time:

Distance	16 cm	22 cm	28 cm	36 cm	45 cm
Position ① (2 tubes)	8500	6800	5000	3500	2500
Recommended duration	35'	45'	1:00	1:25	2:00
Position ② (3 tubes)	14500	10500	7500	5000	3900
Recommended duration	20'	30'	40'	1:00	1:15

Position yourself in front of the appliance and adjust the angle of the Lumino Plus.

- The distance between your face and the screen should be from 15 to 50 cm. Adapt the latter according to your preferences and comfort for your eyes.
- You do not have to sit directly in front of the appliance.
- You can space out the duration of your light therapy session.

Do not look at the lamp directly. During the session, you can read, write, telephone, etc.

Exposure must be regular and during a prolonged period. During the darkest days of the year, repeat the treatment for at least 7 consecutive days, according to individual requirements.

After use, stop the machine by turning the button located on the front of the appliance to position ③ and disconnect it.

If during the session of exposure to light you feel eye pain or other discomfort, immedi-

ately stop using Lumino Plus and consult your doctor.

MAINTENANCE AND STORAGE

Cleaning

Disconnect the appliance and allow it to cool down.

Use a cloth or wet sponge to clean the outside surface of the lamp.

Do not use scouring pads, abrasive products or aggressive detergents such as acetone or petrol to clean the appliance.

Storage

Store the unit in a cool and dry place away from children.

Replacing the tubes

Other than the tubes, the appliance does not contain any other parts to replace. The tubes are designed for approximately 8,000 hours of use (which corresponds to several years of normal operation). The tubes must be replaced with similar tubes.

To replace the tubes, please follow the instructions below:

- Use a screwdriver (not provided) to loosen the 4 screws located on the rear of the product. [■ 5]
- Then remove the reflector cover as shown in the diagram. [■ 6]
- Remove the end clamps and remove the latter as shown in the diagram. [■ 7]
- When you have replaced the tube, put the reflector cover back into place before using the Lumino Plus again.

TECHNICAL CHARACTERISTICS

Voltage: 230V
Frequency: 50Hz
Consumption: 108W

FAQ

⊖ No light

- ⊖ No power
- ⊖ Maybe the appliance is not connected correctly or your power supply is faulty
- ⊖ Master switch off
- ⊖ Switch on the master switch

⊖ One side of the product is no longer operating

- ⊖ Defective neon
- ⊖ Change the neon
- ⊖ Neon at end of service life
- ⊖ Change the neon

⊖ The light continues to blink 20 seconds after being lit

- ⊖ Incorrectly placed neon
- ⊖ Reposition the neon

ADVICE ON WASTE DISPOSAL

The packaging consists entirely of environmentally-friendly materials that can be taken to your local waste sorting centre to be used as secondary materials. Cardboard can be placed in a paper collection container. Packing film must be taken to your local sorting and recycling centre.

When you are no longer using the device, dispose of it in an environmentally responsible manner and in compliance with the statutory requirements.

LIMITED WARRANTY

LANAFORM guarantees that this product will be free from any material and workmanship defects for a period of two years from the date it was purchased, other than as specified below.

The LANAFORM warranty does not cover damage caused by normal wear and tear on this product. In addition, the warranty on this LANAFORM product does not cover damage caused as a result of any abusive, inappropriate or incorrect use, accident, attachment of any unauthorised accessory, modification of the product or any other situation of whatever kind, beyond the control of LANAFORM.

LANAFORM shall not be held liable for any type of incidental, consequential or special damage.

Any implied warranties of the product's fitness for purpose are limited to a period of two years from the original purchase date on presentation of a copy of the proof of purchase.

Upon receipt, LANAFORM will repair or replace your device, as appropriate, and return it to you. The warranty may only be used at the LANAFORM Service Centre. Any maintenance and repair work carried out on this product by any person other than the LANAFORM Service Centre will invalidate this warranty.

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