
LANAFORM

Back Massager

EN ENGLISH	MASSAGE CUSHION	SK SLOVENČINA	MASÁŽNA PODLOŽKA
FR FRANÇAIS	COUVRE-SIÈGE MASSANT	HU MAGYAR	MASSZÍROZÓPÁRNA
NL NEDERLANDS	MASSAGEKUSSEN	SL SLOVENŠČINA	MASAŽNA BLAZINA
DE DEUTSCH	MASSAGESITZ	HR HRVATSKI	MASAŽNI JASTUK
ES ESPAÑOL	COJÍN DE MASAJE	RU РУССКИЙ ЯЗЫК	МАССАЖНАЯ ПОДУШКА
IT ITALIANO	CUSCINO MASSAGGIANTE	BG БЪЛГАРСКИ	МАСАЖНА СЕДАЛКА
PL POLSKI	PODUSZKA MASUJĄCA	SV SVENSKA	MASSAGEDYNA
CS ČEŠTINA	MASÁŽNÍ PODUŠKA		

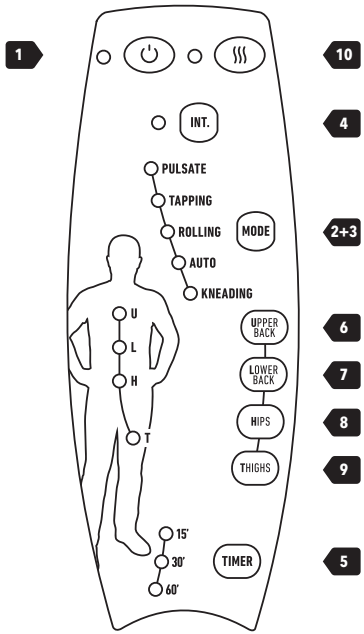☒ **The photographs and other images of the product in this manual and on the packaging are as accurate as possible, but do not claim to be a perfect reproduction of the product.**

Thank you for choosing the Back Massager from LANAFORM, a massage device combining quality and simplicity.

The Back Massager is a massaging seat cover equipped with five motors. These create vibrations which help reduce tension and muscular pain and soothe tired muscles. Comfortable and simple to use, the Back Massager offers a very pleasant vibrating massage. This device is equipped with a function giving you the option to target the massage area and it also has an automatic function which gives you a harmonious sequence of different massage options. In addition, it is operated with a user-friendly control which fits comfortably into the palm of your hand. Back Massager can be used in the office, at home,...

⚠ **Please read all the instructions prior to using your Back Massager, in particular these few basic safety instructions.**

Only use this appliance for the purposes described in this manual.

Before using the device, remember to remove the silica sachet used to combat humidity in the air.

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

Please do not use accessories which are not recommended by LANAFORM or which are not supplied with this appliance.

If the power supply lead is damaged, it must be replaced by a similar lead available from the supplier or its after-sales service.

Do not use this appliance if the plug is damaged, if it is not working correctly, if it has been dropped or if it is damaged or if it has fallen into water. If this is the case, have the appliance inspected and repaired by the supplier or its after-sales department.

Do not carry this appliance by its supply lead or use this lead as a handle.

Always disconnect the appliance from the mains after use or before cleaning it.

Should the appliance fall into water, unplug it immediately before retrieving it.

An electrical appliance should never be left plugged in and unsupervised. Disconnect it from the mains if you are not using it.

Keep the power lead away from hot surfaces.

Do not use this product before going to bed. Massage has a stimulating effect which may inhibit sleep.

Do not use this massage appliance while you are asleep.

Never allow anything to fall into or insert any object (drawing pin or other metal clip) into any of the openings on this device.

Never use this appliance in a room where aerosol products (sprays) are used or in a room where oxygen is administered.

Do not use this appliance under a cover or under a cushion. Excessive heat may cause a fire, electrocution of the individual or other injuries.

When unplugging the appliance, please make sure that the "On/Off" button is switched off before removing the plug.

If you are worried about your health, consult a doctor before using this appliance.

If you feel any pain while using this appliance, stop using it immediately and consult your doctor.

Never use this appliance on a part of the body which is swollen or inflamed or if you have any form of skin rash.

This product is a non-professional massage appliance designed to soothe tired muscles. Do not use this appliance as a replacement for medical care, it is intended exclusively for family use.

How to use

Details

- Massage by vibrations
- 5 motors
- 4 vibration modes plus 1 automatic mode
- 4 massage intensities
- 4 massage areas (upper back, lower back, hips and thighs)
- 1 "heating" function
- Timer (15, 30 or 60min.)

Description of the Back Massager

The strap system

Back Massager has a system of straps enabling you to fix it to any chair. Simply pass the straps around the chair back, bring them together and fix them to each other using the two adhesive ends.

Control device 1

01 ON/OFF

02 4 different massages

03 1 automatic mode

04 4 vibration intensities

05 3 massage durations

06 Upper back massage

07 Lower back massage

08 Hips massage

09 Thighs massage

10 "Heating" function

Using of the Back Massager

ON/OFF

Back Massager operates via an adaptor. First connect the end of the power supply cable to the female input connector located on the right-hand side of the seat. Then plug the adaptor into a power socket. Once connected to the electrical socket, the Back Massager is ready to use.

You can also plug your cushion into your car using the 12 Volt adaptor (provided) and the cigarette lighter. Please remember that the massage cushion must not be used when driving.

Starting or stopping the appliance is controlled by simply pressing the ON/OFF button (diagram, point 1).

Manual mode

Thanks to this mode, you have the option to select:

- One type of massage from the four offered: Mode 1, 2, 3, 4 (diagram, point 2).
- One of the four vibration intensities (diagram, point 4).
- The duration of your massage (15, 30 or 60 minutes) (diagram, point 5).
- The massage area according to your wishes or needs (diagram, points 6, 7, 8 and 9).
- The "heating" function which works to relax your muscles (diagram, point 10).

Ⓞ You may not feel the heat through clothing such as a jumper or jacket.

Automatic mode (diagram, point 3)

This mode consists of a sequence of different massages. This function is valuable since it offers a varied massage. The "TIMER", "INT." and "HEAT" functions and the desired massage area can also be used in automatic mode.

Maintenance

Storage

Put the Back Massager in its box in a dry place.

Avoid it coming into contact with sharp edges and other pointed objects which might cut or pierce the fabric surface.

In order to avoid any breakage of the power lead, do not wind it around the cushion.

Do not hang the cushion up by the manual control lead or the power lead.

Cleaning

Disconnect the cushion from the power socket and wait for a few minutes before cleaning it.

Clean it with a soft, slightly damp sponge.

Never let water or any other liquid come into contact with the cushion.

Do not immerse it in liquid to clean it.

Never use abrasive detergents, brushes or a thinning agent to clean it.

Advice regarding the disposal of waste

All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.

When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law.

Limited warranty

LANAFORM guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM's control.

LANAFORM may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM's Service Centre. Any attempt to maintain this product by a person other than LANAFORM's Service Centre will render this guarantee void.

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PL © Jeśli macie na sobie Państwo sweter lub marynarkę, to możecie nie poczuć ciepła.