LANAFORM®

# Turbo Mass

Vibrating Massager

Vibrating pillow

Coussin vibrant

Trillend kussen

# TURBO MASS - USER MANUAL

## Introduction

Thank you for choosing the Lanaform Turbo Mass, which is practical and easy to use at home and in your car. The vibrations of the Turbo Mass soothe away fatigue and help improve blood circulation. We hope this product will meet with your entire satisfaction.

Please read the following instructions carefully before using the Turbo Mass.

## Important safety instructions

**Read all the instructions before use: incorrect use of this appliance may cause injury or damage.**

- Only connect the Turbo Mass to a socket with the correct nominal voltage for the appliance.
- It is important to check regularly that the cable is undamaged, at each end in particular. If the cable supplied is damaged, please contact your retailer.
- Never use this appliance if the cable or plug is damaged. If it does not function normally, or falls on the floor or into water, contact your retailer.
- Never leave the appliance unsupervised while in use.
- As with all electrical equipment, it is not advisable to use the Turbo Mass during a thunderstorm.
- Keep the cable away from all hot or sharp objects.
- Disconnect the appliance before cleaning it.
- Do not pull on the electrical cable in order to remove the plug.
- Do not dismantle or attempt to repair the appliance yourself, but contact your retailer instead.
- Do not use or store this appliance in a bathroom or any other damp place.

## Caution


- If your appliance smells of burning or the electrical socket smokes, or other problems appear, stop using it. If there is an accident, disconnect the appliance and contact your retailer.
- It is inadvisable to use this equipment on:
  - patients suffering from a tumour;
  - pregnant or menstruating women;
  - persons suffering from cardiac problems or wearing a pacemaker;
  - persons currently undergoing treatment or displaying unusual symptoms.
- The Turbo Mass has not been designed for use by persons (including children) with reduced mental or sensory capacities or with a lack of experience and knowledge, unless they are under the supervision of a person and have been given instructions on using the appliance.
- Keep this appliance out of the reach of children and disabled persons. Keep young children under supervision to ensure they do not play with the appliance.

## Instructions for use

- Connect your appliance. Press the ON/OFF button to switch it on. (vibromassage)
- Press the ON/OFF button once more to switch the appliance off.

- Remove the plug from the socket after use.

### The positions we suggest for use:
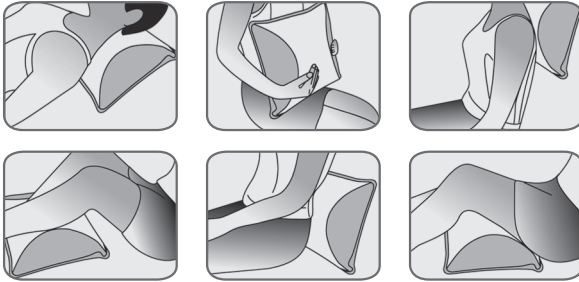
### Maintenance and conservation

- Use a damp rag to clean your appliance and leave it to dry in a well ventilated place.
- Do not use cleaning products on it.
- Do not place your appliance near a fire or where the temperature is very high.
- Do not place your appliance near a source of water.

### Parameters

Voltage: 12 V

Current: 1500 mA

### LIMITED GUARANTEE

LANAFORM® guarantees that this product is free from any defect in material or manufacture for a period of two years from its date of purchase, with the exception of the items detailed below.

The guarantee on this LANAFORM® product does not cover damage caused as a consequence of wrongful use or misuse, accident, the attachment of an unauthorized accessory, a modification of the product, or any other condition of any nature whatsoever outside the control of LANAFORM®.

LANAFORM® shall not be held responsible for any type of accessory, consecutive or special damage.

All implicit guarantees of the suitability of the product are limited to a period of two years from the initial date of purchase.

On receiving your Turbo Mass, LANAFORM® will repair or replace it as appropriate, and will return it to you all expenses paid. Guarantee work is carried out only by the LANAFORM® Service Centre. Any servicing done to this product by any person other than the LANAFORM® Service Centre annuls the present guarantee.

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