LANAFORM®

Stim +

Stimulation Equipment

Rechargeable

Rechargeable

Oplaadbaar

INTRODUCTION TO THE PRODUCT

Thank you for purchasing the **Stim+** by LANAFORM®, the low frequency stimulation system. The Stim + is part of the new generation of low frequency appliances intended for electro-stimulation and muscle building.

PLEASE READ ALL THE INSTRUCTIONS PRIOR TO USING YOUR «STIM+», IN PARTICULAR THESE FEW BASIC SAFETY INSTRUCTIONS:

- Only use this appliance for the purposes described in this manual.
- This appliance is not intended for use by persons (including children) whose physical, sensory or mental capabilities are impaired, or by persons lacking experience or knowledge, unless they have had the benefit of prior supervision and instruction in using the device from a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- Please do not use accessories which are not recommended by LANAFORM® or which are not supplied with this appliance.
- If the power lead is damaged, it must be replaced by a special or similar lead available from the supplier or after-sales department.
- Do not use this appliance if the plug is damaged, if it is not working correctly, if it has been dropped or if it is damaged or if it has fallen into water. If this is the case, have the appliance inspected and repaired by the supplier or its after-sales department.
- Should the appliance fall into water, unplug it immediately before retrieving it.
- Do not transport this appliance by carrying it by its power lead and do not use this lead as a handle.
- Always disconnect the appliance from the mains after use or before cleaning it.
- An electrical appliance should never be left plugged in and unsupervised. Disconnect it from the mains if you are not using it.
- Keep the power lead away from hot surfaces.
- Never let anything fall into and never insert anything into one of the air vents.
- Never use this appliance if the ventilation holes are obstructed.
- Never pour any liquid into the openings of this product.
- Do not use this appliance within 2 hours before or after eating or before going to sleep. As the stimulating effect of the massage could disrupt food digestion or delay the onset of sleep.

- Do not use this massage appliance while you are asleep.
- Never use this appliance in contact with water or in a damp location. This could result in the risk of electric shock.
- This appliance is only suitable for indoor use.
- Never use this appliance in a room where aerosol products (sprays) are used or in a room where oxygen is administered.
- Do not use this appliance under a cover or under a cushion. Excessive heat may cause a fire, electrocution of the individual or other injuries.
- If you are worried about your health or if you're under a medical treatment, consult a doctor before using this appliance.
- If you feel any pain while using this appliance, stop using it immediately and consult your doctor.
- Never use this appliance on a part of the body that is swollen or inflamed or if you have any form of skin rash.
- Do not use this appliance if you are pregnant.
- This product is a non-professional massage device, intended exclusively for home use. Do not use this appliance as a substitute for proper medical care.
- Should a fault occur, please do not attempt to repair the «Stim +» yourself. It must be returned to the supplier or after-sales department.
- Do not use the device when you are driving. This could cause a car accident or other problems due to the stimulation.
- Do not use this device after applying lotion, cream or spray (medical or cosmetic).
- It is not recommendable to use the **Stim+** for more than 20 minutes on the same muscle per session. The number of sessions recommended per day should not exceed two.

SPECIAL HEALTH PRECAUTIONS

- Imperatively consult a doctor before using the appliance in the following cases.
 - Nobody suffering from an acute illness, malignant tumour, infectious disease, cardiac condition;
 - Pregnant women;
 - Nobody with a high fever;

- If you have abnormal blood pressure;
 - If you suffer from abnormal skin or dermatology problems;
 - If you are taking treatment, have suffered an accident or are abnormally sensitive.
 - Nobody who has recently undergone a trauma or surgical operation (within the last 6 months)
- Please do not use this product if you have undergone any surgical organ transplant operation.
 - If your body or skin reacts in an abnormal manner, or if you detect redness or irritation, cease use of the appliance and consult your doctor. Follow the instructions of your doctor.
 - Do not place the electrode on any wound or other physical injury. This could cause infection or aggravate the condition of the injuries.
 - Do not use the appliance near to your heart, above the neck, on the head, mouth or genitals, or near to any areas affected by a skin disorder, etc. Use of the device on these areas could cause an accident or discomfort.

THE BENEFITS OF LOW FREQUENCY TREATMENT:


What is a low frequency exactly?

Frequency refers to the number of electrical stimulations per second for the body and is expressed in Hertz (Hz).

A special current constantly circulates within all living objects. The human body is clearly not an exception to this rule. All of the parts of our bodies, including the heart, brain, our muscles, nerves, etc. emit an imperceptible and very weak current, known as bioelectricity. This current plays a key role in the normal functioning of the body.

The low frequency current causes the muscle to contract and relax. When a muscle is relaxed, large quantities of blood are injected. When a muscle contracts, the metabolised blood is expelled. This action encourages the healthy circulation of blood and therefore contraction and relaxation.

DESCRIPTION OF THE APPLIANCE:

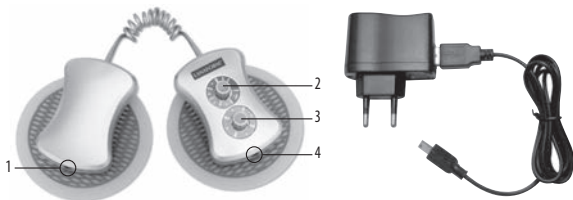1. Charger plug: socket to connect the electrical cable to charge the appliance using the mains supply or USB.
2. Intensity setting knob: turn the switch clockwise to activate the device, the intensity increases on a scale (from 1 to 5) and decreases as you go down the scale, until the device is finally deactivated (OFF position).
3. Mode selection knob: the appliance is programmed with four modes:
 1. Warm-up programme
 2. Beginners programme
 3. Intermediate programme
 4. Advanced programme
4. Connector socket for additional pads (not provided)

METHOD OF USE

First use

1. Prior to the first use, charge the "Stim+" for three hours. The red and green indicator lights light up.
2. When the green indicator lights up, the appliance is fully charged.

Comment: Only charge your battery when completely flat to improve its life cycle.

DIAGRAM OF ELECTRODE POSITIONS

First charge the Stim+.

Warning: Do not forget to remove the protective films from the electrodes prior to the first use.

When using the **Stim+** and when positioning the electrodes, you are recommended to use two electrodes per group of muscles. Electrodes represent positive and negative poles.

We do not recommend placing the two electrodes on either side of your body part. For optimum and safe use, we recommend that you position and connect the electrodes as shown on the diagrams enclosed.

Warning: never apply the electrodes to the head or opposite each other: i.e. do not place the two electrodes on either side of your body part.

GENERAL USE

1. Place the electrodes: Apply the electrodes on the body parts you intend to build up as shown on the diagram.
2. Turn the **Stim+** ON: Turn the intensity knob (knob 2) clockwise from the OFF position. The appliance will start to operate at this point.
3. Select your mode: once you have switched on your **Stim+**, you should select a mode from among the 4 programmes available (warm-up, beginner, intermediate, advanced). Each programme corresponds to the level of the user.
Turn the knob clockwise to select a mode (1 to 4).
4. Set intensity:
Then select the intensity by turning the "Intensity" button clockwise; the intensity increases progressively from 1 to 5.
Turn the intensity knob anti-clockwise; the current will gradually decrease.

We recommend you start with a weak current and increase progressively the first few times you use the appliance. Do not increase the current too rapidly. Once you are used to the stimulation you can increase the intensity progressively.

Always «warm-up» (intensity 1) for a few minutes before each session.

5. Duration of the session: a duration of 20 minutes per session per electrode position is recommended, with one or two sessions per day.

If you wish to modify the position of the adhesive electrode during a session, switch the Intensity knob to the OFF position. You may then reactivate the appliance, to avoid an electric shock.

6. If you wish to stop the session during a treatment, turn the Intensity button to the OFF position, and remove the adhesive electrodes when the appliance is stopped.

MAINTENANCE AND STORAGE

CLEANING

1. Cleaning the main unit:

Use a dry cloth without detergent to remove dirt from the *Stim+*.

2. Cleaning the adhesive sheets

Simply hold the unplugged electrode under clear running water, rinse and rub the adhesive side carefully to remove dust.

3. Storing the adhesive electrodes.

Apply the electrodes to the plastic backing.

STORAGE

- Store the appliance out of the reach of children.
- Do not store the main unit in a location subject to high temperatures, high levels of humidity or direct sunlight. Place it in its box in a dry place.
- Do not dismantle the unit without authorisation; do not drop or knock the main unit.

TROUBLESHOOTING

Breakdown	Checks	Solutions
Low stimulation or the user feels prickling	Check that the electrode is in contact with the skin	Position the adhesive disk of the electrode firmly on the skin
	Check that the surface of the adhesive disk is not dirty	Clean or replace the adhesive disk
	Check that no grease is present on your skin	Clean your skin
The appliance stops its cycle before the end	Check that the battery is charged	Charge the Stim+
The skin treated becomes red, itches or prickles	Check that you are not allergic to the adhesive disk	Check your allergy record for your skin. If your allergic reaction is minor, simply change the position of the adhesive disk or shorten the duration of the treatment. If your allergic reaction is severe, cease all use of the appliance or use anti-allergy medication
Even at maximum intensity, the stimulation is too weak or the user lacks any sensation	Is your skin dirty? Are the electrodes positioned correctly?	Re-use the appliance after cleaning your skin with a wet flannel
The electrodes do not stick anymore		Dip them in water for a few seconds

ADVICE REGARDING THE DISPOSAL OF WASTE

- All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.
- When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law.

LIMITED WARRANTY

LANAFORM® guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM® guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM® product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM®'s control.

LANAFORM® may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM® will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM®'s Service Centre®. Any attempt to maintain this product by a person other than LANAFORM®'s Service Centre will render this guarantee void.

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